BIOTEST

Your reliable home test kits series

Pregnancy Test Kit (One-Step)

Product Code 9 555223 208272
For in vitro diagnostic and self testing use.
Store at 2-30ºC

INTENDED USE

BIOTEST One-Step Pregnancy Test Kit is a rapid chromatographic immunoassay for the qualitative detection of human chorionic gonadotropin (hCG) in urine to aid in the early detection of pregnancy.


CONTENTS

- Test Device
- Desiccant
- Instruction for Use

PRECAUTIONS & WARNINGS

1. Do not use the test after the expiration date.
2. If the package is not completely sealed do not use the test.
3. The test should be performed in a well-lit area.
4. Use the test device immediately after opening it.
5. Do not touch the result window. This could affect results and may also impose personal hazards.
6. The pouch contains a Silica Gel pack to absorb humidity. Do not open the pack. Throw it away with the remaining of the test.
7. At the end of the test, wrap every thing you have used in a plastic bag and throw away in the bin. Do not forget to wash your hands properly.


TEST PROCEDURE

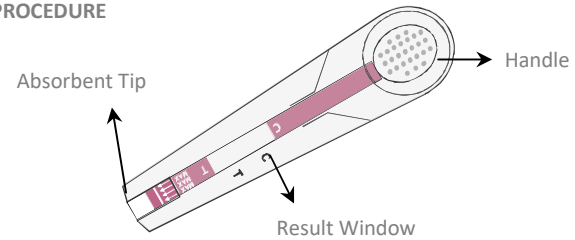


1. Read instructions carefully before use. It is preferred to perform the test on the first morning urine because the first morning urine contains the highest amount of hCG. Check that the pouch is completely sealed then open the pouch.

2. Hold the test kit at its holder and place the other tip under the urine stream for 5 seconds.
3. Place the test kit on a flat surface and wait for the result.
4. Read the result in the result window between 1- 5minutes. Test result should not be read after 10 minutes.


HOW TO READ THE RESULTS

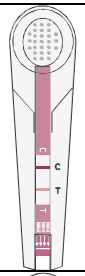

PREGNANT
Positive results : two bands.
The presence of two red lines in the result window means you are probably pregnant. One line may be lighter than the other; they do not have to match. Check with your doctor to confirm pregnancy.

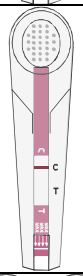

NOT PREGNANT
Negative results : one band in the (C) area.
The presence of one red line in the control **area** (C) means you are probably not pregnant.

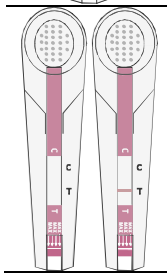

INVALID RESULT
Invalid results : no bands or one band in the (T) area.
The result is invalid if no red line appears in the control area (C), even if a line appears in the test area (T). You should repeat the test with a new test.

WHAT THE RESULTS MEAN

A positive result indicates that your urine contains hCG (the pregnancy hormone) at a concentration more than or equal to 25mIU/ml which is our device sensitivity and you can assume you are pregnant. You should consult your doctor who is best able to further advise you on proper prenatal care. If you are taking medication, inform your doctor immediately and make sure it is safe for your baby.

A negative result means that no hCG has been detected and you can assume you are not pregnant. If a week passes and your menstrual period is still delayed, you should do another test. If the result is still negative in the later test, there is little chance that you are pregnant. However, because there could be health-related reasons why you are not menstruating, you should see your doctor without delay.

If your menstrual cycle is irregular, it can be difficult to know the day you have missed your period. If your test indicated you are not pregnant and your menstrual cycle does not begin within a few days, you may have miscalculated your cycle and should use a second test after two days.

LIMITATION OF THE TEST

The test is not reusable and works only if the instructions are followed precisely. Although the test is highly accurate in detecting pregnancy, a low incidence of false results can occur. Check with your doctor if you get unexpected or inconsistent test results. Certain health conditions can also cause a false or irregular test result.

QUESTIONS & ANSWERS

Q. How does the test work?
A. When a woman is pregnant, her body produces a special pregnancy hormone known as hCG (human Chorionic Gonadotropin). This hormone appears in the urine as early as the first day of a missed menstrual cycle. A positive result indicates that your urine contains hCG and you can assume you are pregnant. A negative result means that no hCG has been detected and you can assume you are not pregnant.

Q. What if the two red lines appear, but one is much lighter than the other, am I pregnant?
A. Yes. You should assume you are pregnant.

Q. Can the result change after a certain length of times?
A. A positive test result might not change after the test is completed and several hours thereafter. However, a negative test result should not be read after 10 minutes.

Q. What should I do if the test result is negative and I still suspect that I am pregnant?
A. You may have miscalculated when your period was due. If your period does not start within a week, repeat the test. If you still get a negative result and your period has not started, you should see your doctor.

Q. What if I am uncertain about the other pink line in the "T" area?
A. It is recommended that you retest three days later using another tests.

IFU-BT-HCG-SH-v1.2018

Manufactured by Atlas Medical
William James House, Cowley Road, Cambridge, CB4 0WX, UK
BIOTEST is a trademark of Bio-gro Sdn Bhd.
www.bio-test.info